




Elemental Scientific Semiconductor Catalog for PerkinElmer Elan

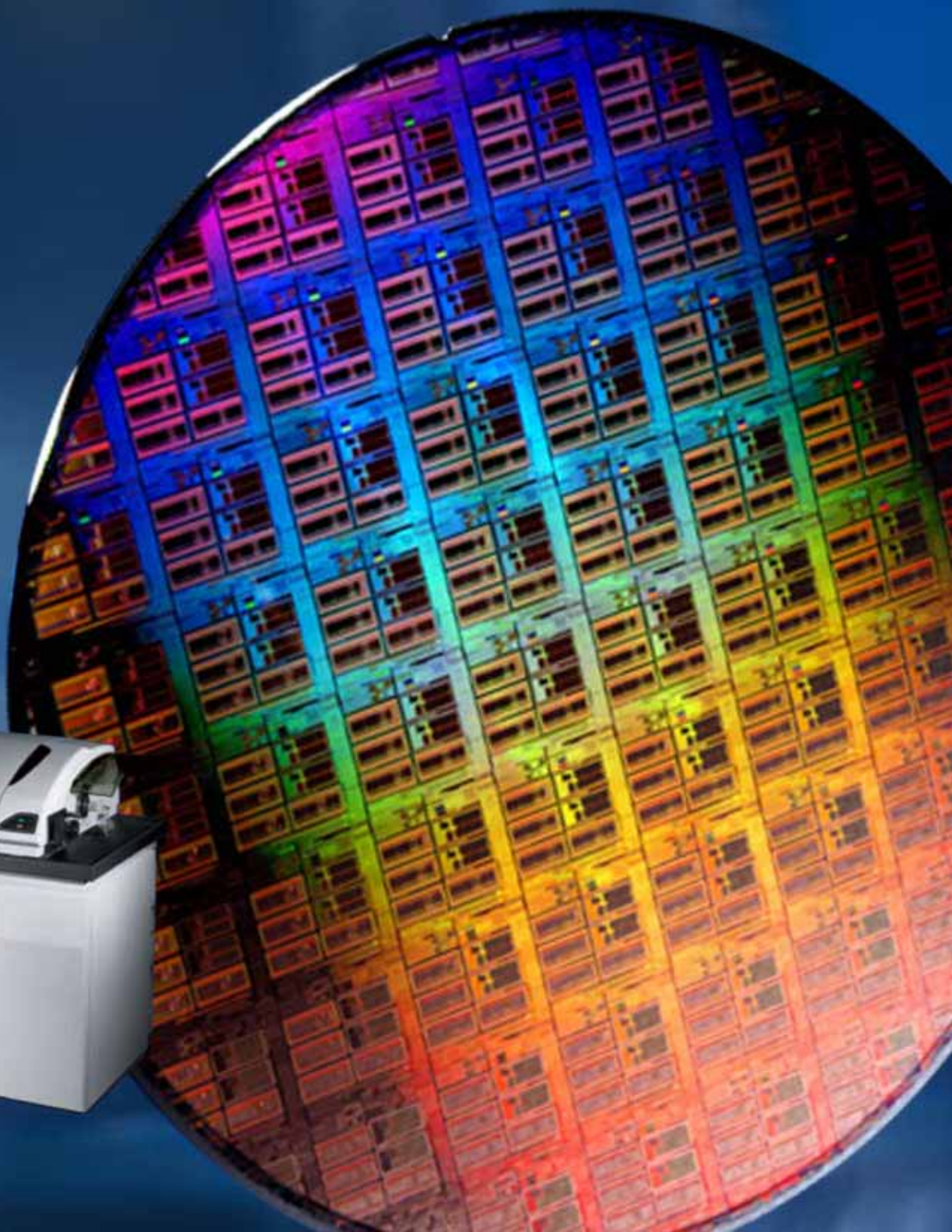

Ordering Information

- You may contact us via phone, fax, email, or through our website.

EPOND S.A.

- CP 389
- CH-1800 Vevey
- Switzerland

Phone:

- +41 21 921 29 41

Fax:

- +41 21 921 44 08

E-Mail:

- info@e-pond.biz

Website:

- <http://www.e-pond.biz>


Table of Contents:

▪ PFA-Pt Sample Introduction Kit .4	
▪ MicroFlow PFA Nebulizers	5
▪ Teflon Sample Probes.....	7
▪ Apex High Sensitivity Upgrade ..	8
▪ Apex Shelves.....	9
▪ PC ³ Peltier Cooled Cyclonic Chamber.....	10
▪ Cyclonic Spray Chambers	10
▪ SC-2 Autosampler.....	11
▪ Autosampler Clean Enclosure .	12
▪ Autosampler Racks.....	13
▪ PFA Vials and Tubes	15
▪ Platinum ICP Cones	16
▪ PFA Pipette Tips	16

High Efficiency PFA MicroFlow Nebulizers

PFA Nebulizers with Integrated Sampling Probes

These nebulizers allow self-aspiration at low sample consumption rates with any autosampler while providing the same benefits of the PFA MicroFlow nebulizer. The integrated sampling probe is constructed by encapsulating a rigid support material between layers of Teflon to provide a chemically resistant probe that is mechanically stable and resilient.

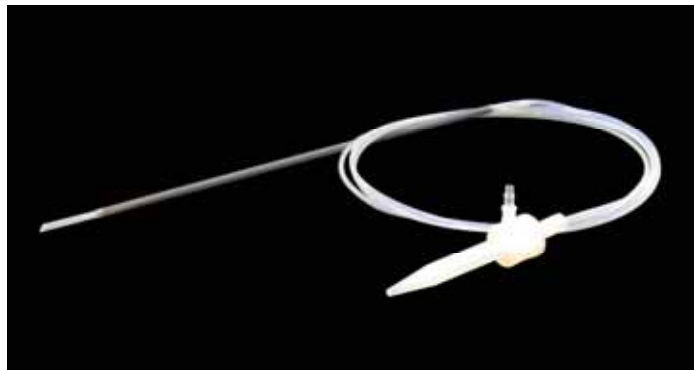

LLL = distance from nebulizer to the probe support (cm). Standard lengths are 80 cm (-080), 100 cm (-100), 120 cm (-120), 150 cm (-150). Probes can be custom-made to your specifications of length and materials at no extra charge.

Carbon Fiber support

Model	ESI SC-series	PE AS-9x	ASX-100	ASX-500	Manual probes
PFA-20	ES-2020-3503-LLL		ES-2020-0103-LLL		ES-2020-0203-LLL
PFA-50	ES-2000-3503-LLL	ES-2000-9503-LLL	ES-2000-0103-LLL	ES-2000-0503-LLL	ES-2000-0203-LLL
PFA-75	ES-2001-3503-LLL	ES-2001-9503-LLL	ES-2001-0103-LLL	ES-2001-0503-LLL	ES-2001-0203-LLL
PFA-100	ES-2002-3503-LLL	ES-2002-9503-LLL	ES-2002-0103-LLL	ES-2002-0503-LLL	ES-2002-0203-LLL
PFA-200	ES-2003-3503-LLL	ES-2003-9503-LLL		ES-2003-0503-LLL	ES-2003-0203-LLL
PFA-400	ES-2005-3503-LLL	ES-2005-9503-LLL		ES-2005-0503-LLL	ES-2005-0203-LLL

PEEK support

Model	ESI SC-series	PE AS-9x	ASX-100	ASX-500	Manual probes
PFA-20	ES-2020-3500-LLL		ES-2020-0100-LLL		ES-2020-0200-LLL
PFA-50	ES-2000-3500-LLL	ES-2000-9500-LLL	ES-2000-0100-LLL	ES-2000-0500-LLL	ES-2000-0200-LLL
PFA-75	ES-2001-3500-LLL	ES-2001-9500-LLL	ES-2001-0100-LLL	ES-2001-0500-LLL	ES-2001-0200-LLL
PFA-100	ES-2002-3500-LLL	ES-2002-9500-LLL	ES-2002-0100-LLL	ES-2002-0500-LLL	ES-2002-0200-LLL
PFA-200	ES-2003-3500-LLL	ES-2003-9503-LLL		ES-2003-0500-LLL	ES-2003-0200-LLL
PFA-400	ES-2005-3500-LLL	ES-2005-9503-LLL		ES-2005-0500-LLL	ES-2005-0200-LLL

Ultem support

Model	ESI SC-series	PE AS-9x	ASX-100	ASX-500	Manual probes
PFA-20	ES-2020-3505-LLL		ES-2020-0105-LLL		ES-2020-0205-LLL
PFA-50	ES-2000-3505-LLL	ES-2000-9505-LLL	ES-2000-0105-LLL	ES-2000-0505-LLL	ES-2000-0205-LLL
PFA-75	ES-2001-3505-LLL	ES-2001-9505-LLL	ES-2001-0105-LLL	ES-2001-0505-LLL	ES-2001-0205-LLL
PFA-100	ES-2002-3505-LLL	ES-2002-9505-LLL	ES-2002-0105-LLL	ES-2002-0505-LLL	ES-2002-0205-LLL
PFA-200	ES-2003-3505-LLL	ES-2003-9503-LLL		ES-2003-0505-LLL	ES-2003-0205-LLL
PFA-400	ES-2005-3505-LLL	ES-2005-9503-LLL		ES-2005-0505-LLL	ES-2005-0205-LLL


Autosampler Probes for ST Nebulizers

ST Nebulizer probes connect directly to any ST-type nebulizer. The 1/4-28 threaded fitting provides a secure, zero-dead-volume connection with no additional fittings required.

LLL = distance from flanged fitting to the probe support (cm).

- 80 cm (-080)
- 100 cm (-100)
- 120 cm (-120)
- 150 cm (-150)

Probes can be custom-made to your specifications of length and materials at no extra charge.

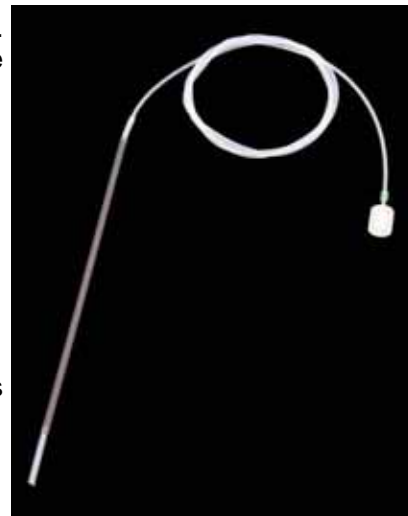

Carbon Fiber Support

Probe I.D.	ESI-SC	PE AS-9x	ASX-100	ASX-500	Manual probes
0.15mm	ES-5037-3150-LLL	ES-5037-9150-LLL	ES-5036-0150-LLL	ES-5037-0150-LLL	ES-5036-2150-LLL
0.20mm	ES-5037-3200-LLL	ES-5037-9200-LLL	ES-5036-0200-LLL	ES-5037-0200-LLL	ES-5036-2200-LLL
0.25mm	ES-5037-3250-LLL	ES-5037-9250-LLL	ES-5036-0250-LLL	ES-5037-0250-LLL	ES-5036-2250-LLL
0.30mm	ES-5037-3300-LLL	ES-5037-9300-LLL	ES-5036-0300-LLL	ES-5037-0300-LLL	ES-5036-2300-LLL
0.50mm	ES-5037-3500-LLL	ES-5037-9500-LLL		ES-5037-0500-LLL	ES-5036-2500-LLL
0.80mm	ES-5037-3750-LLL	ES-5037-9750-LLL		ES-5037-0750-LLL	ES-5036-2750-LLL

PEEK Support

Probe I.D.	ESI-SC	PE AS-9x	ASX-100	ASX-500	Manual probes
0.15mm	ES-5017-3150-LLL	ES-5017-9150-LLL	ES-5036-0150-LLL	ES-5017-0150-LLL	ES-5016-2150-LLL
0.20mm	ES-5017-3200-LLL	ES-5017-9200-LLL	ES-5036-0200-LLL	ES-5017-0200-LLL	ES-5016-2200-LLL
0.25mm	ES-5017-3250-LLL	ES-5017-9250-LLL	ES-5036-0250-LLL	ES-5017-0250-LLL	ES-5016-2250-LLL
0.30mm	ES-5017-3300-LLL	ES-5017-9300-LLL	ES-5036-0300-LLL	ES-5017-0300-LLL	ES-5016-2300-LLL
0.50mm	ES-5017-3500-LLL	ES-5017-9500-LLL		ES-5017-0500-LLL	ES-5016-2500-LLL
0.80mm	ES-5017-3750-LLL	ES-5017-9750-LLL		ES-5017-0750-LLL	ES-5016-2750-LLL

Ultem Support

Probe I.D.	ESI-SC	PE AS-9x	ASX-100	ASX-500	Manual probes
0.15mm	ES-5057-3150-LLL	ES-5057-9150-LLL	ES-5056-0150-LLL	ES-5057-0150-LLL	ES-5056-2150-LLL
0.20mm	ES-5057-3200-LLL	ES-5057-9200-LLL	ES-5056-0200-LLL	ES-5057-0200-LLL	ES-5056-2200-LLL
0.25mm	ES-5057-3250-LLL	ES-5057-9250-LLL	ES-5056-0250-LLL	ES-5057-0250-LLL	ES-5056-2250-LLL
0.30mm	ES-5057-3300-LLL	ES-5057-9300-LLL	ES-5056-0300-LLL	ES-5057-0300-LLL	ES-5056-2300-LLL
0.50mm	ES-5057-3500-LLL	ES-5057-9500-LLL		ES-5057-0500-LLL	ES-5056-2500-LLL
0.80mm	ES-5057-3750-LLL	ES-5057-9750-LLL		ES-5057-0750-LLL	ES-5056-2750-LLL

Apex High Sensitivity Inlet System

Fully-integrated, high-sensitivity sample introduction systems for ICP spectroscopy. Designed to increase signal intensity while reducing relative background contamination for ICPMS analyses. When combined with a self-aspirating MicroFlow PFA nebulizer even aggressive and concentrated chemicals may be analyzed at very high sensitivity with low blank contamination.

Apex Q for PE Elan

High-sensitivity, low-flow sample introduction system for ICPMS. Improves sample transport efficiency while maintaining low blanks. Recommended to be used with the Apex Shelf. Includes PFA-100 and PFA-ST nebulizers, integrated peristaltic drain pump, and connectors to adapt to the PE ELAN.

Apex Q for ELAN: **ES-4354**



Apex HF for PE Elan

All PFA sample flow path. High-sensitivity, low-flow HF-resistant sample introduction system for ICPMS. Recommended to be used with the Apex shelf. Includes PFA-100 and PFA-ST nebulizers, integrated peristaltic drain pump, and connectors to adapt to the PE ELAN.

Apex HF for ELAN: **ES-4454**



Spiro Teflon Membrane Desolvator

Designed for use with or without the Apex system. Innovative design allows for direct analysis of a variety of samples, including high concentration Si samples and volatile organic solvents. The Spiro is an extremely stable membrane desolvator compatible with many types of high dissolved solids samples. The Spiro removes Si on-line as SiF_6 after addition of HF to the sample, allowing the determination of trace metal impurities in Si matrices without ICPMS cone clogging or off-line sample preparation. The Spiro combined with the Apex system is useful for direct analysis of neat organic chemicals, such as methanol, acetone and other solvents. Since the Spiro removes volatile species from the sample aerosol stream, it is not recommended for the determination of volatile analytes, including boron.

Spiro Membrane Desolvator: **4599-3000**



The Spiro attaches directly to the Apex.




Apex Shelves for Elan ICPMS

Polypropolyne shelves for direct mounting of the Apex system to the Elan ICPMS injector. A notch in the bottom of each shelf fits directly over the ridge in the built-in Elan tray.

Apex Shelf for PE Elan with bayonet torch

for PE Elan 6100, DRC,
9000 before mid-2005: **ES-2999-5431**

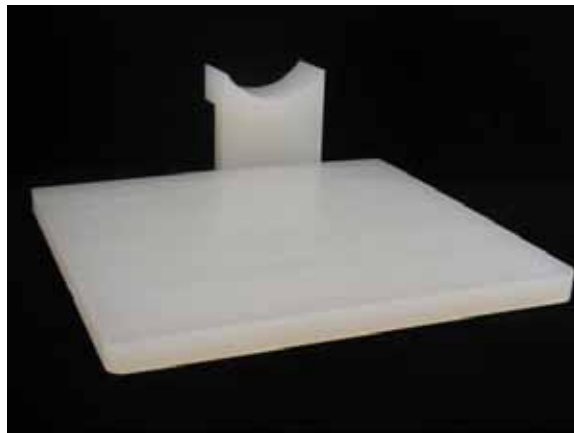

Apex Shelf for PE Elan with cassette torch

for PE Elan newer
than mid-2005: **ES-2999-5421**



PC³ Peltier Cooler for Cyclonic Spray Chambers

The PC³ is a small and robust Peltier cooler for ICPMS cyclonic spray chambers. The cyclonic spray chamber offers the advantages of low-memory effects, fast rinse-out, and higher sample transport efficiency. The PC³ cools the outer walls of the cyclonic spray chamber to a constant temperature enhancing long-term signal stability, reducing polyatomic ion interferences, such as oxides, and reducing solvent loading when volatile organic solvents are analyzed.

PC³ kit with Quartz spray chamber and addition gas port: **ES-4152-1111**

PC³ kit with Polypro spray chamber and addition gas port: **ES-4152-2111**



The PC3 can be adapted to a wide variety of ICP injectors, including ESI PFA injector base and quartz ball joint injector. The PC3 shown here is connected to the ELAN injector normally used with the PE Ryton spray chamber.



O-ring-free Cyclonic Spray Chambers

Cyclonic spray chambers offer the advantages of low-memory effects, fast rinse-out, and higher sample transport efficiency. O-ring free connections to the nebulizer and injector eliminate a possible source of contamination. These spray chambers can be used with any 6mm nebulizer, including MicroFlow PFA Nebulizers, PolyPro Nebulizers, and glass concentric nebulizers.

Quartz spray chamber with addition gas port: **ES-4152-1111**

Polypro spray chamber with addition gas port: **ES-4152-2111**



A one piece ferrule nut holds any 6mm nebulizer in place

*PC³ cyclonic spray chambers may be used alone or with the PC³ Peltier Cooler.

SC-2 Ultraclean Optically Encoded Autosampler

Automatic sampling while maintaining sample integrity and purity. Improved control of self-aspiration, reducing "vapor-locks" and useful for low flow applications and organic solvents analysis. Demonstrably improves rinse-out and background over any other autosampler available for the PE Elan. Recently improved for chemical resistance to high acid concentrations in enclosed environments. Fully compatible with ELAN and Optima software.

SC-2 Autosampler
with IEEE converter: **SC-8102C**



SC-FAST high purity sample injection

The SC-FAST integrates fully with the SC-2 autosampler to ensure the fastest sample throughput of any commercial sample introduction system. Through use of an injection valve system, most of the uptake and rinse time can be eliminated, increasing throughput up to 2—3 fold without loss of analytical performance.

SC-FAST upgrade: **SC-0500-03**




SC-FAST can be integrated to the bottom of the SC-2 for PE Elan ICPMS.

Semiconductor rinse stations and delivery

Dual pressurized containers deliver the rinse solution to the SC Autosampler dual-flowing rinse stations without the use of a peristaltic pump, eliminating a source of contamination. All wetted parts are made from PFA Teflon to ensure the cleanest possible rinse.

PFA Rinse Delivery System: **SC-0304**

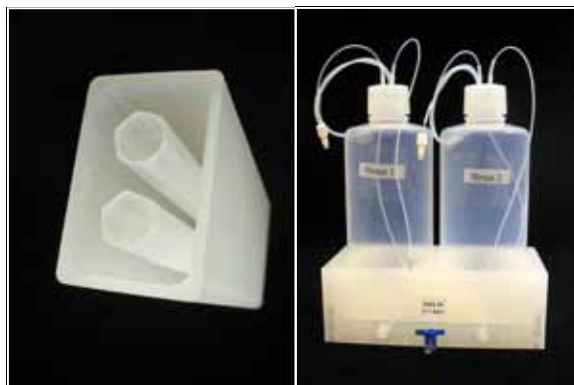

Upgrade options for SC-2 Autosampler

Sample enclosure for SC-2 autosampler

Controls evaporation and environmental contamination. Only the sample area is enclosed allowing free access to the rest of the instrument. The large flexible door opens easily for access to the sample racks and locks securely closed.

Sample enclosure: **SC-1207**

Please contact ESI if you are interested in HEPA / ULPA filtered sample environments for the SC Autosampler



SC-2 STAND

The combination of the SC-2 autosampler with STAND option provides an optimized sample changer for the Elan ICPMS. The stand positions the SC-2 at an ideal height next to the Elan ICPMS, minimizing the distance between the sample and the nebulizer for faster uptake and rinse-out. The SC-2 STAND also eliminates the need for a separate autosampler table and includes two shelves with storage space.

SC-2 STAND: **SC-1210**



The STAND option can be used either with the SC-2 autosampler along or combined with the SC-FAST sample injection system.

Sample Racks for SC-2




MR-21-14 21 position 3x7 mini rack. For "GFAA" style vials.




MR-00-01 Cover for Micro Racks.



MT-24G 24 position 4x6 GeMeTec rack. Compatible with GeMeTec automated VPD system.

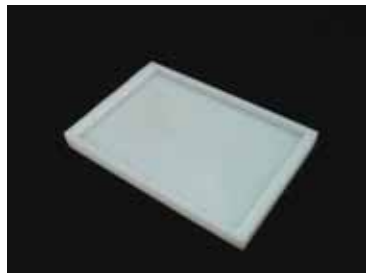

MT-24G-AP GeMeTec adapter plate. Adapts GeMeTec rack to microtiter plate locations on the SC-2 autosampler.

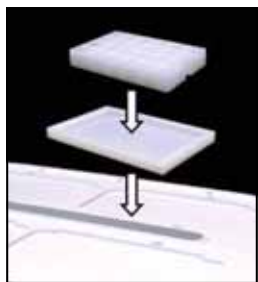

The SC autosampler readily accommodates a wide variety of sample racks for both small and large volume sample analysis. Many other sizes and types of racks are available. Please contact ESI for racks compatible with vials that you prefer to use.

Sample Racks for SC-2




ST10-28 10 position standards rack. For 25 mL PFA bottles.



SR2-21-27 21 position 3x7 super rack. For 25 ml PFA bottles



SR2-12-50 12 position 2x6 super rack. For 125 ml PFA bottles

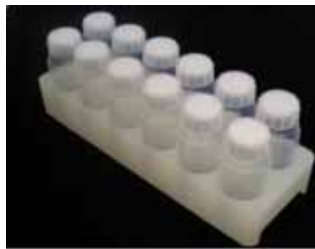

The SC autosampler readily accommodates a wide variety of sample racks for both small and large volume sample analysis. Many other sizes and types of racks are available. Please contact ESI for racks compatible with vials that you prefer to use.

Sample Vials and Bottles for SC-2

V-14-0311C 1.5mL PFA "GFAA" style vial
(press-on cap not included)



V-14-0312C 2mL PFA "GFAA" style vial
(press-on cap not included)



V-14-0314C 4mL PFA "GFAA" style vial,
(press-on cap not included)



V-14-0310 PFA push on cap for "GFAA" style vials.



V-27-0360 25mL PFA bottle with PFA seal and PFA screw cap



V-50-0360 125mL high purity PFA bottle with seal and screw cap.
Recommended for DI water and high purity semiconductor-grade chemicals due to favorable sample volume to surface area ratio.

ICPMS Cones for PE Elan

Platinum Cones

ES-3000-1711 Pt Sampler cone
 ES-3000-1712 Pt Skimmer cone

Nickel Cones

ES-3000-1709 Ni Sampler cone
 ES-3000-1710 Ni Skimmer cone



Elan 6000/6100/9000/DRC Cones

PFA Pipette Tips for Eppendorf type Pipettors

100 μ l PFA Pipette Tips

Quantity of 10, pre-cleaned: **ES-7000-0101**
 Quantity of 96, pre-cleaned: **ES-7000-0111**



1000 μ l PFA Pipette Tips

Quantity of 10, pre-cleaned: **ES-7000-1001**
 Quantity of 60, pre-cleaned: **ES-7000-1011**

